


# CT/MR Contour Prostate Template Sets

SETM1204-02-XX, SETM1204-04-XX, SETM1206-02-XX, SETM1206-04-XX

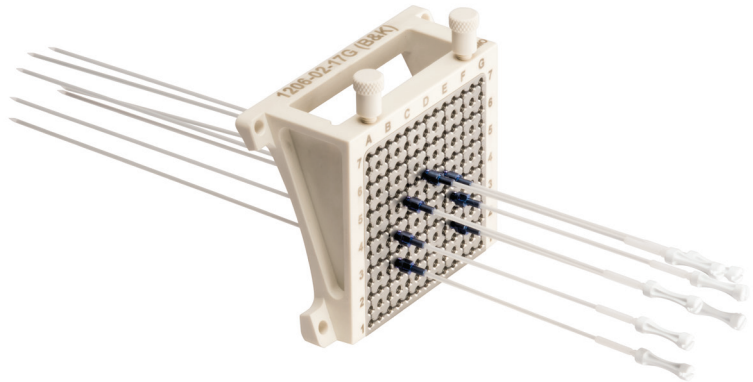

Facilitates interstitial needles for the treatment of prostate cancer.

## Prostate

### Simplicity

The CT/MR Contour Prostate Template facilitates the placement of interstitial needles for the treatment of prostate cancer. All templates are laser engraved with letters to denote the different vertical columns of needle channels and with numbers to denote the different horizontal rows of needle channels. The numbering/lettering layout corresponds to the type of ultrasound system being used.

### Flexibility

The CT/MR Contour Prostate Template is available to accommodate needles of different gauges, which are securely held in place by color coded needle collets associated with the appropriate needle size. Both small and large template sizes are provided. The CT/MR Contour Prostate Template is curved to fit the sloping anatomy of the perineum. There are four suture holes provided for optional fixation of the template to the perineum.


### Reliability

The CT/MR Contour Prostate Template is constructed of PPSU (Polyphenylsulfone) and the template inserts and needle collets are constructed of titanium for continuous use and steam sterilization.

### Features

- Contoured to fit the sloping perineum
- Expected lifetime is 5 years
- MR conditionality allows for 1.5T or 3T
- Individual locking collet design enhances discrete needle targeting

- Available for variety of needle sizes: 15.5G, 17G, 18G, 1.65 mm, 2 mm
- Provided with four suture holes allowing the template to be directly sutured to the perineum
- Resistance-free bore holes for smooth needle entry
- Durable construction

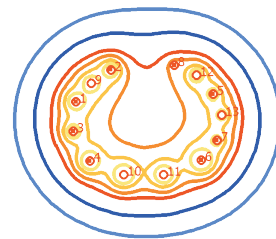

Isodose



Template needle collets



Eckert & Ziegler

Mick Radio-Nuclear  
Instruments, Inc.

An Eckert & Ziegler BEBIG Company

# CT/MR Contour Prostate Template Sets

## SETM1204-02-XX (Small Matrix), for CIVCO Steppers, BK labeling scheme\*

Qty.	Item	Part number
1	CT/MR Contour Prostate Template PPSU Body, small matrix	
1	CT/MR Contour Prostate Template Titanium Insert, small matrix	
1	Titanium Collet Wrench	LCZ1205-02
1	Mounting Pegs, (2 pcs)	LCZ1205-03
1	Mounting Bracket Lock Knobs, (2 pcs)	LCZ1205-12
<i>Choose one CT/MR Contour Prostate Template, small matrix, BK based on needle size (XX)</i>		
1	CT/MR Contour Prostate Template, small matrix, BK, 15.5G	LCZ1204-02-15.5
1	CT/MR Contour Prostate Template, small matrix, BK, 17G	LCZ1204-02-17
1	CT/MR Contour Prostate Template, small matrix, BK, 17G, Eckert & Ziegler BEBIG	LCZ1204-02-EZ-17
1	CT/MR Contour Prostate Template, small matrix, BK, 18G	LCZ1204-02-18
1	CT/MR Contour Prostate Template, small matrix, BK, 1.65 mm	LCZ1204-02-1.65
1	CT/MR Contour Prostate Template, small matrix, BK, 1.65 mm, Eckert & Ziegler BEBIG	LCZ1204-02-EZ-1.65
1	CT/MR Contour Prostate Template, small matrix, BK, 2.0 mm	LCZ1204-02-2.0
<i>Choose one Template Needle Collets Set (36 pcs) based on needle size (XX)</i>		
6	Template Needle Collets (Titanium), gray, 15.5G, (6 pcs)	LCZ1105-01-15.5
6	Template Needle Collets (Titanium), blue, 17G, (6 pcs)	LCZ1105-01-17
6	Template Needle Collets (Titanium), seafoam, 17G, Eckert & Ziegler BEBIG, (6 pcs)	LCZ1105-01-EZ-17
6	Template Needle Collets (Titanium), gold, 18G, (6 pcs)	LCZ1105-01-18
6	Template Needle Collets (Titanium), magenta, 1.65 mm, (6 pcs)	LCZ1105-01-1.65
6	Template Needle Collets (Titanium), bronze, 1.65 mm, Eckert & Ziegler BEBIG, (6 pcs)	LCZ1105-01-EZ-1.65
6	Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)	LCZ1105-01-2.0

## SETM1204-04-XX (Small Matrix), for CIVCO Steppers, Siemens labeling scheme\*

Qty.	Item	Part number
1	CT/MR Contour Prostate Template PPSU Body, small matrix	
1	CT/MR Contour Prostate Template Titanium Insert, small matrix	
1	Titanium Collet Wrench	LCZ1205-02
1	Mounting Pegs, (2 pcs)	LCZ1205-03
1	Mounting Bracket Lock Knobs, (2 pcs)	LCZ1205-12
<i>Choose one CT/MR Contour Prostate Template, small matrix, Siemens based on needle size (XX)</i>		
1	CT/MR Contour Prostate Template, small matrix, Siemens, 15.5G	LCZ1204-04-15.5
1	CT/MR Contour Prostate Template, small matrix, Siemens, 17G	LCZ1204-04-17
1	CT/MR Contour Prostate Template, small matrix, Siemens, 17G, Eckert & Ziegler BEBIG	LCZ1204-04-EZ-17
1	CT/MR Contour Prostate Template, small matrix, Siemens, 18G	LCZ1204-04-18
1	CT/MR Contour Prostate Template, small matrix, Siemens, 1.65 mm	LCZ1204-04-1.65
1	CT/MR Contour Prostate Template, small matrix, Siemens, 1.65 mm, Eckert & Ziegler BEBIG	LCZ1204-04-EZ-1.65
1	CT/MR Contour Prostate Template, small matrix, Siemens, 2.0 mm	LCZ1204-04-2.0
<i>Choose one Template Needle Collets Set (36 pcs) based on needle size (XX)</i>		
6	Template Needle Collets (Titanium), gray, 15.5G, (6 pcs)	LCZ1105-01-15.5
6	Template Needle Collets (Titanium), blue, 17G, (6 pcs)	LCZ1105-01-17
6	Template Needle Collets (Titanium), seafoam, 17G, Eckert & Ziegler BEBIG, (6 pcs)	LCZ1105-01-EZ-17
6	Template Needle Collets (Titanium), gold, 18G, (6 pcs)	LCZ1105-01-18
6	Template Needle Collets (Titanium), magenta, 1.65 mm, (6 pcs)	LCZ1105-01-1.65
6	Template Needle Collets (Titanium), bronze, 1.65 mm, Eckert & Ziegler BEBIG, (6 pcs)	LCZ1105-01-EZ-1.65
6	Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)	LCZ1105-01-2.0



Eckert & Ziegler BEBIG SagiNova® afterloader

# CT/MR Contour Prostate Template Sets

## SETM1206-02-XX (Large Matrix), for CIVCO Steppers, BK labeling scheme\*

Qty.	Item	Part number
1	CT/MR Contour Prostate Template PPSU Body, large matrix	
1	CT/MR Contour Prostate Template Titanium Insert, large matrix	
1	Titanium Collet Wrench	LCZ1205-02
1	Mounting Pegs, (2 pcs)	LCZ1205-03
1	Mounting Bracket Lock Knobs, (2 pcs)	LCZ1205-12
<i>Choose one CT/MR Contour Prostate Template, large matrix, BK based on needle size (XX)</i>		
1	CT/MR Contour Prostate Template, large matrix, BK, 15.5G	LCZ1206-02-15.5
1	CT/MR Contour Prostate Template, large matrix, BK, 17G	LCZ1206-02-17
1	CT/MR Contour Prostate Template, large matrix, BK, 17G, Eckert & Ziegler BEBIG	LCZ1206-02-EZ-17
1	CT/MR Contour Prostate Template, large matrix, BK, 18G	LCZ1206-02-18
1	CT/MR Contour Prostate Template, large matrix, BK, 1.65 mm	LCZ1206-02-1.65
1	CT/MR Contour Prostate Template, large matrix, BK, 1.65 mm, Eckert & Ziegler BEBIG	LCZ1206-02-EZ-1.65
1	CT/MR Contour Prostate Template, large matrix, BK, 2.0 mm	LCZ1206-02-2.0
<i>Choose one Template Needle Collets Set (36 pcs) based on needle size (XX)</i>		
6	Template Needle Collets (Titanium), gray, 15.5G, (6 pcs)	LCZ1105-01-15.5
6	Template Needle Collets (Titanium), blue, 17G, (6 pcs)	LCZ1105-01-17
6	Template Needle Collets (Titanium), seafoam, 17G, Eckert & Ziegler BEBIG, (6 pcs)	LCZ1105-01-EZ-17
6	Template Needle Collets (Titanium), gold, 18G, (6 pcs)	LCZ1105-01-18
6	Template Needle Collets (Titanium), magenta, 1.65 mm, (6 pcs)	LCZ1105-01-1.65
6	Template Needle Collets (Titanium), bronze, 1.65 mm, Eckert & Ziegler BEBIG, (6 pcs)	LCZ1105-01-EZ-1.65
6	Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)	LCZ1105-01-2.0

### Required Items

**SETM1204-02-EZ-17, SETM1204-04-EZ-17, SETM1206-02-EZ-17, SETM1206-04-EZ-17**

Qty.	Item	Part number
1	Disposable Steel Needle Set	SET0052
20	Transfer Tube, D = 3 mm, L = 1000 mm	LAF1000

### Required Items

**SETM1204-02-EZ-1.65, SETM1204-04-EZ-1.65, SETM1206-02-EZ-1.65, SETM1206-04-EZ-1.65**

Qty.	Item	Part number
1	Plastic Needle Set	SET0012
20	Easy Click Transfer Tube, D = 3 mm, L = 1000 mm	LAG1000

## SETM1206-04-XX (Large Matrix), for CIVCO Steppers, Siemens labeling scheme\*

Qty.	Item	Part number
1	CT/MR Contour Prostate Template PPSU Body, large matrix	
1	CT/MR Contour Prostate Template Titanium Insert, large matrix	
1	Titanium Collet Wrench	LCZ1205-02
1	Mounting Pegs, (2 pcs)	LCZ1205-03
1	Mounting Bracket Lock Knobs, (2 pcs)	LCZ1205-12
<i>Choose one CT/MR Contour Prostate Template, large matrix, Siemens based on needle size (XX)</i>		
1	CT/MR Contour Prostate Template, large matrix, Siemens, 15.5G	LCZ1206-04-15.5
1	CT/MR Contour Prostate Template, large matrix, Siemens, 17G	LCZ1206-04-17
1	CT/MR Contour Prostate Template, large matrix, Siemens, 17G, Eckert & Ziegler BEBIG	LCZ1206-04-EZ-17
1	CT/MR Contour Prostate Template, large matrix, Siemens, 18G	LCZ1206-04-18
1	CT/MR Contour Prostate Template, large matrix, Siemens, 1.65 mm	LCZ1206-04-1.65
1	CT/MR Contour Prostate Template, large matrix, Siemens, 1.65 mm, Eckert & Ziegler BEBIG	LCZ1206-04-EZ-1.65
1	CT/MR Contour Prostate Template, large matrix, Siemens, 2.0 mm	LCZ1206-04-2.0
<i>Choose one Template Needle Collets Set (36 pcs) based on needle size (XX)</i>		
6	Template Needle Collets (Titanium), gray, 15.5G, (6 pcs)	LCZ1105-01-15.5
6	Template Needle Collets (Titanium), blue, 17G, (6 pcs)	LCZ1105-01-17
6	Template Needle Collets (Titanium), seafoam, 17G, Eckert & Ziegler BEBIG, (6 pcs)	LCZ1105-01-EZ-17
6	Template Needle Collets (Titanium), gold, 18G, (6 pcs)	LCZ1105-01-18
6	Template Needle Collets (Titanium), magenta, 1.65 mm, (6 pcs)	LCZ1105-01-1.65
6	Template Needle Collets (Titanium), bronze, 1.65 mm, Eckert & Ziegler BEBIG, (6 pcs)	LCZ1105-01-EZ-1.65
6	Template Needle Collets (Titanium), green, 2.0 mm, (6 pcs)	LCZ1105-01-2.0

### Required Items

**SETM1204-02-2.0, SETM1204-04-2.0, SETM1206-02-2.0, SETM1206-04-2.0**

Qty.	Item	Part number
20	Plastic Needle, sterile, D = 2.0 mm, L = 250 mm, 15G	LLA250-K15
5	CT/MR Obturator for 15G Plastic Needle, sterile, L = 250 mm	LLH04-250B15
20	Easy Click Transfer Tube 6F, D = 3 mm, L = 1000 mm	LAG1000-6F

### Optional Items

Qty.	Item	Part number
1	Basket with lid for cleaning and sterilizing of small parts, 121 x 73 x 36 mm**	LLZ07-01
1	Sterilization Box, 381 x 254 x 38 mm***	LLZ05-01

\* The CT/MR Contour Prostate Template has been tested for compatibility with the following accessories and external devices based on stated nominal dimensions. Users should validate that their particular needle is compatible with the template and collet prior to use:  
- Plastic & metal needle diameters that can be used with the corresponding template include: 15.5G, 17G, 18G, 1.65 mm (16G) and 2.0 mm  
- CIVCO Steppers with BK and Siemens labeling scheme  
CIVCO is a trade name/trademark of CIVCO Medical Instruments Co. Inc. BK Ultrasound is a trade name/trademark of Analogic Corporation or its affiliated companies. Siemens is a trade name/trademark of Siemens & Halske Aktiengesellschaft.  
\*\* Manufactured by Kögel GmbH  
\*\*\* Manufactured by PST Plastic Sterilizing Trays Corporation  
Mick and SagiNova are registered trademarks of Eckert & Ziegler BEBIG GmbH and its subsidiaries.  
The mentioned products are not available in all markets. Please contact your local Eckert & Ziegler BEBIG representative for more information.

**Manufacturer:**

**Mick Radio-Nuclear Instruments, Inc.  
An Eckert & Ziegler BEBIG Company**  
521 Homestead Avenue  
Mount Vernon, NY 10550  
USA

Phone +1 914 667 3999  
Fax +1 914 665 8834  
sales@micknuclear.com

**Regional Sales, Marketing and Service:**

**Europe, Middle East, Africa,  
Latin America, Asia Pacific**

**Eckert & Ziegler  
BEBIG GmbH**  
Robert-Rössle-Str. 10  
13125 Berlin  
Germany

Phone +49 30 94 10 84 130  
Fax +49 30 94 10 84 112  
info@bebig.com

**North America**

**Mick Radio-Nuclear Instruments, Inc.  
An Eckert & Ziegler BEBIG Company**  
521 Homestead Avenue  
Mount Vernon, NY 10550  
USA

Phone +1 914 667 3999  
Fax +1 914 665 8834  
sales@micknuclear.com

www.bebig.com  
www.micknuclear.com