


# H.A.M. Applicators (Flaps)

LFO01-03 - LFO01-24

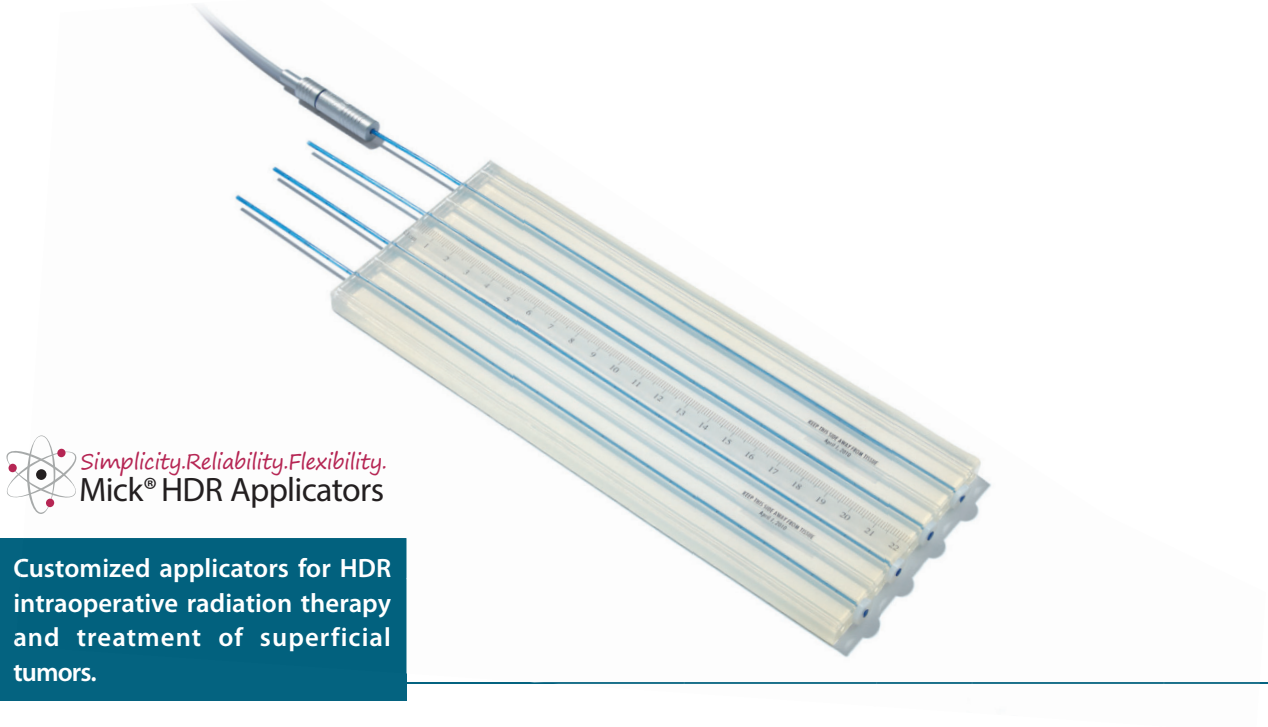

 *Simplicity. Reliability. Flexibility.*  
Mick® HDR Applicators

Customized applicators for HDR intraoperative radiation therapy and treatment of superficial tumors.

## Skin, Mucosa, Intraoperative

### Various Treatment Sites

The Harrison-Anderson-Mick Applicators are used for the treatment of skin surface, mucosa and tumor bed/organ surface (intraoperative).

### Flexibility

The H.A.M. Applicator consists of a flexible pad of silicone rubber that is 8 mm thick. An array of catheters are embedded parallel to each other spaced 10 mm apart while a consistent source-to-tissue distance of 5 mm is maintained. Because of its flexibility, the H.A.M. Applicator easily conforms to the shape of surfaces to which it is applied.

### Customized Solution

H.A.M. Applicators are available in different sizes. Multiple numbers of channels from 3 up to 24 provide the desired treatment area coverage. The H.A.M. Applicators are for single use only.

### Features

- Various treatment options
- Guarantee uniform dose distribution
- Customized sizes
- Each H.A.M. Applicator can be ordered also separately
- Single use

### LFO01-03 - LFO01-24

Quantity	Item	Part number
1	H.A.M. Applicator, 3-channel, 6 Fr, without embedded catheters	LFO01-03
1	H.A.M. Applicator, 4-channel, 6 Fr, without embedded catheters	LFO01-04
1	H.A.M. Applicator, 6-channel, 6 Fr, without embedded catheters	LFO01-06
1	H.A.M. Applicator, 9-channel, 6 Fr, without embedded catheters	LFO01-09
1	H.A.M. Applicator, 12-channel, 6 Fr, without embedded catheters	LFO01-12
1	H.A.M. Applicator, 15-channel, 6 Fr, without embedded catheters	LFO01-15
1	H.A.M. Applicator, 18-channel, 6 Fr, without embedded catheters	LFO01-18
1	H.A.M. Applicator, 21-channel, 6 Fr, without embedded catheters	LFO01-21
1	H.A.M. Applicator, 24-channel, 6 Fr, without embedded catheters	LFO01-24


### Required Items

Quantity	Item	Part number
1	Easy Click Transfer Tube, D = 3 mm, L = 1000 mm	LAG1000
20	Flexible Catheter Single Leader, D = 1.65 mm, L = 300 mm	LLA300-F
1	Cutter for Flexible Catheters <sup>7</sup>	LAH06-01

### Optional Items

Quantity	Item	Part number
1	Set of X-Ray Markers, L = 325 mm, Code = 01-20	SET0034
*****	Fixation Balls, d = 1.6 mm, D = 7 mm, (10 pcs)	LLZ06-05
1	Opening Forceps for Fixation Balls <sup>2</sup>	LAH05-02
***** Depending on the numbers of channels to be used		

### Isodose

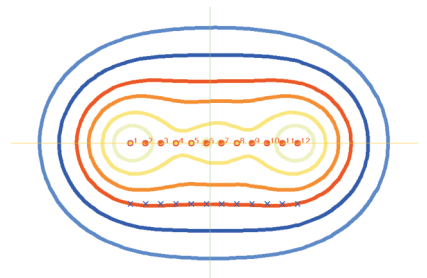


### Easy Click Transfer Tube



### Cutter for Flexible Catheter



### X-Ray Markers

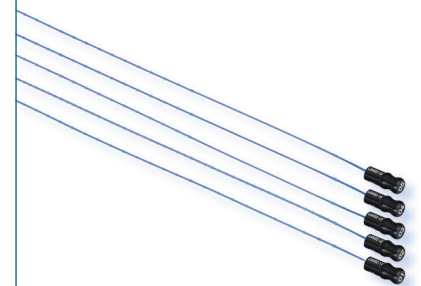

<sup>2</sup> Manufactured by Otto Leibinger GmbH

<sup>7</sup> Manufactured by Pressmaster GmbH

